

REMOVABLE BUCKY INSERTS

Product Advantages

- Converts existing bucky cassette trays to fit any cassette size wireless panel
- Lightweight design with a handle and Twin-Lock provides superior panel protection
- Can be used in both the table and upright wall stand
- Can have customized cutouts for battery, LEDs and power button



Replacing a table and wall bucky is expensive. The material and labor costs to upgrade the equipment can add up fast. The **Bucky Insert** is an adaptor that updates the cassette tray at a low cost without the need for a field service engineer. The outside dimensions of the **Bucky Insert** mimic the outside dimensions of the previous oversized DR panel, while the internal dimensions of the **Bucky Insert** match that of an ANSI standard wireless 14x17in. detector. It can be customized with cutouts for the battery, LEDs and power button. This way the panel can always remain safe and secure in this protective case. The single ergonomic handle provides a sure confident grip, and decreases the chance for accidental drops when moving from the table to the wall. Milled from a single piece of HDPE, the holder is strong, easy to clean, and features our Twin-Lock dual locking system.


Product Details




Tube Side



Rear Side

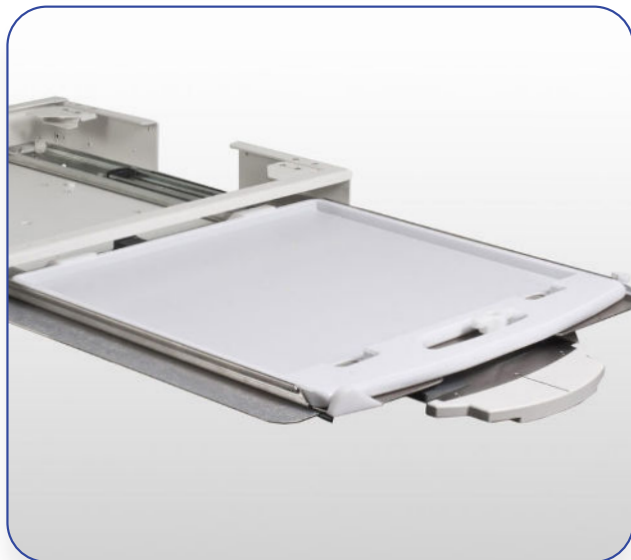

Bucky Insert in table equipment

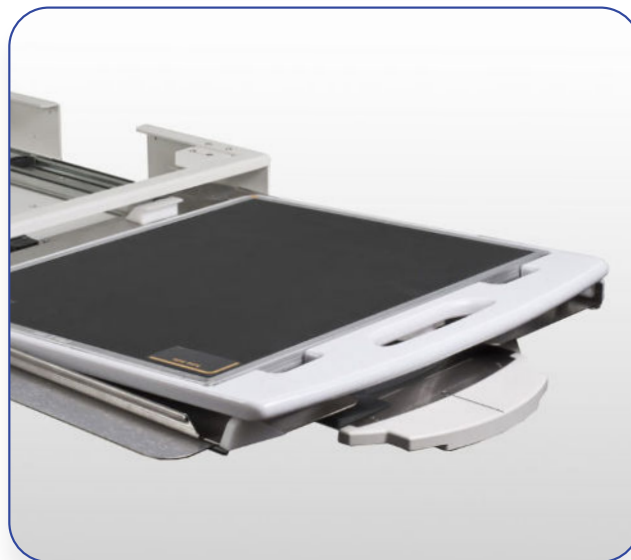

Bucky Insert with DR panel in table equipment

Technical Data

BI50G

Product Dimensions

19.25x18.875x.875in. (48.895x47.943x2.22cm)

Weight

1.6 lbs. (.73 kg)

BI55G

Product Dimensions

18.875x18.937x.875in. (47.943x48.099x2.22cm)

Weight

2.15 lbs. (.98 kg)



Reinalmaging CE

Safer • Protective • Preferred Radiology Accessories



www.reinaimaging.com



customer_service@reinaimaging.com